

3211857

<http://www.synergyhealthplc.com>

Certificate of Irradiation

Date Issued: 16-Nov-2015

FR01S11506960-1-1

This is to certify that Synergy Health Marseille has where appropriate delivered an irradiation process in accordance with:

EN ISO 11137-1:2015 Sterilisation of Health Care Products

EN ISO 9001:2008 Quality Management System

EN ISO 13485:2012 Quality System - Medical Devices

ECP ENTEGRIS CLEANING PROCESS

Le Millénaire

395 rue Louis Lépine

34000 MONTPELLIER

FRANCE

DIVERSEY

clearklaus Bi-Spore RTU VH26S
Batch N° - ENT26381 15 310
DIVERSEY, 18.12.2015

Order Information

Account Number:	101999
Synergy Health Sales Part Reference:	1023395
Customer Reference Number:	CF151420 - 12/11/2015
Product Description:	224514P - BOTTLES 100 ML + STOPPERS
Validation Reference:	S11215019
Quantity Received:	1
Customer Minimum Specification kGy:	25.0
Customer Maximum Specification kGy:	40.0

Irradiation Data

Date and Time of Irradiation:	15-NOV-2015 01:21
Calculated Minimum Dose kGy:	26.1
Calculated Maximum Dose kGy:	33.6


Irradiation Release Authorised By Synergy Health plc

Processing Site: M.I.N. 712 - Arnavaux, , Marseille Cedex 14, 13323 Phone No: +33 (0) 4 91 214 214

Registered Office: M.I.N. 712 - Arnavaux, 13323 Marseille Cedex 14, FRANCE

N° TVA: FR59 343 092 540

Page 1 of 1

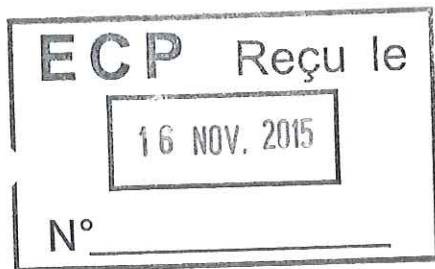

3211179

<http://www.synergyhealthplc.com>

Certificate of Irradiation

Date Issued: 16-Nov-2015

FR01S11506958-1-1

This is to certify that Synergy Health Marseille has where appropriate delivered an irradiation process in accordance with:

EN ISO 11137-1:2015 Sterilisation of Health Care Products

EN ISO 9001:2008 Quality Management System

EN ISO 13485:2012 Quality System - Medical Devices

ECP ENTEGRIS CLEANING PROCESS

Le Millénaire

395 rue Louis Lépine

34000 MONTPELLIER

FRANCE

DIVERSEY

Charkleus Bi-Spore RTU VH26S

Batch N° ENT86381 15 310

26 NOV 84, 18.12.2015

Order Information

Account Number:	101999
Synergy Health Sales Part Reference:	1023397
Customer Reference Number:	CF151419 - 12/11/2015
Product Description:	224516P - KITS BIDONS 5L
Validation Reference:	179848
Quantity Received:	1
Customer Minimum Specification kGy:	15.0
Customer Maximum Specification kGy:	30.0

Irradiation Data

Date and Time of Irradiation:	15-NOV-2015 13:14
Calculated Minimum Dose kGy:	20.1
Calculated Maximum Dose kGy:	26.8

Irradiation Release Authorised By Synergy Health plc

Processing Site: M.I.N. 712 - Arnavaux, , Marseille Cedex 14, 13323 Phone No: +33 (0) 4 91 214 214

Registered Office: M.I.N. 712 - Arnavaux, 13323 Marseille Cedex 14, FRANCE

N° TVA: FR59 343 092 540

Page 1 of 1

Customer DIVERSEY LTD	Management Quality Manual n° MMQ0000 (contractual)
---------------------------------	--

Product : Clearklens Bi-Spore RTU VH26S	Control report n° : CERT 15666
Specification reference : NA	Purchase order n° : 4701742670
Analyst : MJA	Date of reception : 11/23/2015

PHYSICAL AND CHEMICAL PROPERTIES CONTROL

Customer's batch n° : ENT26381 15 310	SYNERGY HEALTH certificate(s) n° : S11506900101
Date of manufacturing : 11/2015	Produced quantity : 160 x 5L
Expiry date : 11/2017	FIDT N° : 1816 ind Z
ECP's batch n° : 26381	Date of analysis : 11/23/2015

Results :				
Solution	Analyzed characteristics	Method of analysis	State in the production	Specification
			Beginning	
Bi-Spore Activator Solution	Appearance (20°C)	NA	Clear colourless liquid	Clear colourless liquid
	pH (20°C)	NA	11,53	9,0 - 12,5
	Specific gravity (20°C)	IDT 0301	1,005	1,004 - 1,020
Bi-Spore Base Solution RTU	Appearance (20°C)	NA	Clear colourless liquid	Clear colourless liquid
	pH (20°C)	NA	2,76	2,5 - 5,0
	Specific gravity (20°C)	IDT 0301	1,000	1,000 - 1,010
Bi-Spore Activator Solution + Base Solution RTU	Appearance (20°C)	NA	Clear yellow liquid	Clear yellow liquid
	pH (20°C)	NA	2,97	2,5 - 5,0
	Specific gravity (20°C)	IDT 0301	1,000	1,000 - 1,010
	Dosage of the materiel chlorine dioxide	IDT 0261	121,5	100 - 180 ppm
Microbiological contagion of the WFI		IDT 0236	4	< 10 cfu/100mL

C : Conform NC : Not Conform
NA : Not Applicable

Conclusion :

☒ CONFORM ☐ NOT CONFORM

Batch released by D.GHERSY, 18.12.2015

Date : 11/30/2015	Edited : E.BENOIT	Checked : D.GHERSY
-------------------	-------------------	--------------------

18 DEC. 2015

Rapport N° : 15/OI/STE/470

Report N° :

Demandeur : ECP

Enquirer: 395, rue Louis Lépine
34000 MONTPELLIER

Rapport d'essai - Essai de stérilité par la méthode de filtration sur membrane
Test report - Membrane filtration sterility testing

selon la Pharmacopée Européenne 8^e édition chapitre 2.6.1
according to the 8th European Pharmacopoeia Edition, chapter 2.6.1

ECHANTILLON(S)/SAMPLE(S)

Désignation : Clearklens Bi-spore Numéro de commande : CF151500
Product name : Order number :
Référence client : 7516197 POE : -
Customer reference : SIP : -
Numéro de lot : ENT26381 15 310 Matériau(x) : -
Batch number : Material : -
Date de réception : mercredi 2 décembre 2015 Donnée de stérilisation : -
Receipt date : Sterilization data :
Date de réalisation : jeudi 3 décembre 2015 Commentaire(s) : 10 bidons poolés
Testing date : Comment(s) :
Nombre d'échantillon(s) : 2
Sample quantity :

PROTOCOLE/PROTOCOL

Volume d'échantillon : 500 ml 3 rinçages de 100 ml de DNP Thio
Sample volume :
Solution neutralisante : DNP + Thiosulfate Volume de rinçage : 3 x 100 ml
Neutralizing solution : Rinsing volume :
Nombre de milieux testés : 2 Volume d'immersion de la membrane : 100 ml
Number of Environment tested : Immersion volume of membranes :

Conditions Conditions	Milieux de culture Environment of culture	Température d'incubation Incubation temperature	Durée d'incubation Incubation period
Bactéries aérobies, levures, champignons Aerobic and fungal	Bouillon Trypcase Soja Tryptone soy solution	22,5 ± 2,5°C	14 jours 14 days
Bactéries anaérobies et aérobies Anaerobic and aerobic	Bouillon Thioglycolate avec Résazurine Thioglycolate Resazurine	32,5 ± 2,5°C	14 jours 14 days

Validation de la méthode
Method validation

09/OI/VAL.STE/017

RESULTATS / RESULTS

Conditions / Milieux de culture Conditions / Media of culture	Examen de la croissance microbienne du milieu Examination of the microbial growth in the media			
	Après 7 jours After 7 days		Après 14 jours After 14 days	
Bactéries aérobies, levures et moisissures / Bouillon Trypcase Soja Aerobic and fungal / Tryptone soy solution	-	Positifs Positive	0	Positifs Positive
	-	Négatifs Negative	1	Négatifs Negative
Bactéries anaérobies et aérobies / Bouillon Thioglycolate avec Résazurine Anaerobic and aerobic / Thioglycolate Resazurine solution	-	Positifs Positive	0	Positifs Positive
	-	Négatifs Negative	1	Négatifs Negative

CONTRÔLES / CONTROLS

Avant (Before)		Après (After)	
Contrôle plan de travail (UFC) Work plan control (CFU)	0	0	
Air sous flux laminaire (UFC) Laminar flow air control (CFU)	0		
		Contrôle de gant (UFC) Glove control (CFU) :	1
		Contrôle stérilité milieux : Media sterility control :	Conforme

CONCLUSION

Les échantillons testés ne présentent pas de croissance microbienne après 14 jours d'incubation
Tested products doesn't shown any microbial development after 14 days of incubation.

Aucun produit ne s'est révélé positif lors de cet essai.

No Product was positive during the test.

Rédigé par : Justine DORAND
Written by : Technicienne Biologiste
Biologist technician

Approuvé par : Audrey BERTHOME
Approved by : Responsable d'essais microbiologiques
Microbiological tests Manager

Date : vendredi 18 décembre 2015

Le présent rapport ne peut être reproduit totalement ou partiellement qu'avec l'approbation du laboratoire et ne concerne que les objets soumis à l'essai.
This report may not be reproduced in whole or part without the approval of the laboratory and only involves the tested products.